

ABC COMBO TEST

FLU A+B & COVID-19 AG

NASAL SWAB

REF RCV-RD0301-01

MDA Reg. No.: IVD5790924-181335

For professional *in vitro* diagnostic use only.**SCAN ME**Please scan
the QR code
for more
information**EASY**
nasal swab
procedure**FAST**
results in
20 minutes**RELIABLE**
test
performance

PRODUCT SPECIFICATION

**Sensitivity**Influenza A: 91.25%
Influenza B: 92.05%
COVID-19: 96.05%**Test time**

20 minutes

**Sample type**Nasal swab
specimen**Specificity**

100%

* vs FDA Emergency Use
Authorized real-time
RT-PCR**Storage condition**

2-30°C

**Shelf life**


24 months

IVD


RESZON® ABC COMBO TEST - FLU A+B & COVID-19 Ag



Gently tear open the pouch and remove the test device

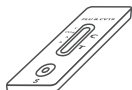

Lay the device on a clean, flat work surface. Label the test device with sample name

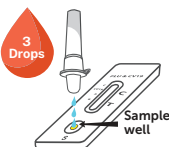

Add 3 drops (~75µl) of the sample solution to the sample well. Make sure there is no air bubble. Sample will start wicking up the membrane.




Read results within 10-20 minutes.

Disregard the results after 30 minutes.
A positive result may be visible at 3 minutes.

Complete reaction time of 20 minutes is required to confirm a negative result.

RESULT INTERPRETATION

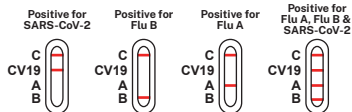

Positive



Negative



Invalid

REAGENTS AND MATERIALS SUPPLIED

- RESZON® ABC Combo Test - Flu A+B & COVID-19 Ag test device (sealed in foil pouch)
- Extraction buffer tube and dropper
- Disposable swab (in sterile packing)
- Biohazard waste bag
- Instruction for Use (product insert)

ORDERING INFORMATION

RESZON® ABC COMBO TEST - FLU A+B & COVID-19 Ag

Product Code	Description	Packing Size
RCV-RD0301-01	RESZON® ABC Combo Test - Flu A+B & COVID-19 Ag (Nasal Swab)	1 test / kit
RCV-RD0301-10	RESZON® ABC Combo Test - Flu A+B & COVID-19 Ag (Nasal Swab)	10 tests / kit
RCV-RD0301-25	RESZON® ABC Combo Test - Flu A+B & COVID-19 Ag (Nasal Swab)	25 tests / kit

Contact us today to enquire: +6012 815 8689 | info@reszonics.com

Reszon Diagnostics International Sdn. Bhd. 201001033978 (917901-W)

A member of Hextar Healthcare Berhad

©2026 Reszon Diagnostics International Sdn. Bhd. All rights reserved. RESZON® is the registered trademark of Reszon Diagnostics International Sdn. Bhd. Doc No: MM-01-04-13-EN; Rev. No.: 04; Effective Date: 2026.08.11
The content of this material is intended for Healthcare Professionals only, not for general public.

www.reszonics.com

Reszon
Diagnostics
INTERNATIONAL

Reszon Diagnostics

