



FILARIASIS IgG4 RAPID TEST KIT

SCAN ME



Please scan
the QR code
for more
information

REF RLF-RD0101-25

MDA Reg. No.: IVDC3603823-154220

For professional *in vitro* diagnostic use only.



Detect specific IgG4 antibodies to

- *Wuchereria bancrofti*
- *Brugia malayi*
- *Brugia timori*



HASSLE-FREE

no special
equipment
needed



EASY

sample preparation
protocol and visual
result interpretation



RELIABLE

test
performance

PRODUCT SPECIFICATION



Sensitivity
96.53%



Test time
10-15 minutes



Sample volume
20 µl



Sample type
Human Whole Blood /
Serum / Plasma



Specificity
99.57%



Storage
2-30°C




Kit
25 tests



Shelf life
24 months

IVD

RESZON® FILARIASIS IgG4 RAPID TEST KIT

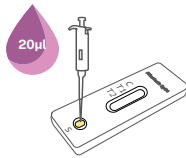

Pipette 20 µl of serum, plasma or whole blood into the sample well.



Add 2 drops of buffer to the same sample well.




Sample will start wicking up the membrane




Read the test result within 10-15 minutes

- 1 Positive result may appear as early as 5-10 minutes.
- Negative result must be confirmed only at 15 minutes.
- Do not interpret the result after 20 minutes.


RESULT INTERPRETATION

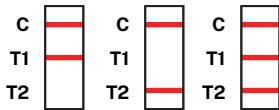

Positive

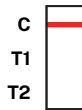

Negative

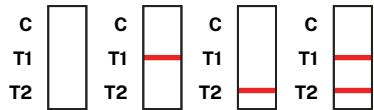

Invalid

REAGENTS AND MATERIALS SUPPLIED

- RESZON® Filariasis IgG4 Rapid Test cassette (individually sealed aluminium pouch with a desiccant)
- Assay buffer
- Instruction for Use (product insert)

ORDERING INFORMATION

RESZON® FILARIASIS IgG4 RAPID TEST KIT

Product Code	Description	Packing Size
RLF-RD0101-25	RESZON® Filariasis IgG4 Rapid Test Kit	25 tests / kit

Contact us today to enquire: +6012 815 8689 | info@reszonics.com

Reszon Diagnostics International Sdn. Bhd. 201001033978 (917901-W)

A member of Hextar Healthcare Berhad

©2026 Reszon Diagnostics International Sdn. Bhd. All rights reserved. RESZON® is the registered trademark of Reszon Diagnostics International Sdn. Bhd. Doc No: MM-01-04-04-EN; Rev. No.: 05; Effective Date: 2026.08.10
The content of this material is intended for Healthcare Professionals only, not for general public.

www.reszonics.com

Reszon
Diagnostics
INTERNATIONAL

Reszon Diagnostics

