

HBsAg RAPID TEST

REF RHB-RD0101-25, RHB-RD0101A-25, RHB-RD0101B-25

MDA Reg. No.: IVDD4340025-214346

For professional *in vitro* diagnostic use only.

SCAN ME



Please scan
the QR code
for more
information



Qualitative detection of Hepatitis B surface antigen (HBsAg) in human whole blood, serum or plasma specimens



HASSLE-FREE

no special
equipment
needed



FAST

results in
20 minutes



EASY

sample preparation
protocol and visual
result interpretation



RELIABLE

test
performance

PRODUCT SPECIFICATION



Sensitivity
99.3%



Test time
15-20 minutes



Sample volume
75 µl



Sample type
Human Whole Blood /
Serum / Plasma



Specificity
99.8%



Storage
2-30°C





Kit
25 tests



Shelf life
24 months

RESZON® HBsAg RAPID TEST

FOR WHOLE BLOOD SAMPLE

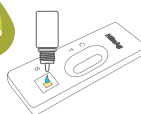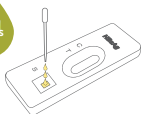

Add 3 drops (approximately 75 µl) of whole blood into the sample well marked "S" at the bottom.

Add 1 drop of buffer to the same sample well.

The results may be read from 15 to 20 minutes

1 Allow the reaction to proceed for 15 minutes, read the results.
Do not read results after 20 minutes.

FOR SERUM AND PLASMA SAMPLE




Add 3 drops (approximately 75 µl) of serum or plasma into the sample well marked "S" at the bottom.

The results may be read from 15 to 20 minutes

1 Allow the reaction to proceed for 15 minutes, read the results.
Do not read results after 20 minutes.


RESULT INTERPRETATION

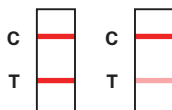

Positive



Negative

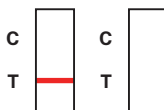

Invalid

REAGENTS AND MATERIALS SUPPLIED

- RESZON® HBsAg Rapid test cassette (individually sealed aluminium pouch with a desiccant)
- Assay buffer
- Optional: Disposable dropper, alcohol swab and safety lancet
- Instruction for Use (product insert)

ORDERING INFORMATION

RESZON® HBsAg RAPID TEST

Product Code	Description	Packing Size
RHB-RD0101-25	RESZON® HBsAg Rapid Test	25 tests / kit
RHB-RD0101A-25	RESZON® HBsAg Rapid Test (with disposable dropper)	25 tests / kit
RHB-RD0101B-25	RESZON® HBsAg Rapid Test (with disposable dropper, alcohol swab and safety lancet)	25 tests / kit

Contact us today to enquire: +6012 815 8689 | info@reszonics.com

Reszon Diagnostics International Sdn. Bhd. 201001033978 (917901-W)

A member of Hextar Healthcare Berhad

©2026 Reszon Diagnostics International Sdn. Bhd. All rights reserved. RESZON® is the registered trademark of Reszon Diagnostics International Sdn. Bhd. Doc No: MM-01-04-36-EN; Rev. No.: 03; Effective Date: 2026.08.10
The content of this material is intended for Healthcare Professionals only, not for general public.

www.reszonics.com

Reszon
Diagnostics
INTERNATIONAL

Reszon Diagnostics

