



HCV RAPID TEST

REF RHC-RD0101-25, RHC-RD0101A-25, RHC-RD0101B-25

MDA Reg. No.: IVDD10681325-211023

For professional *in vitro* diagnostic use only.

SCAN ME



Please scan
the QR code
for more
information



Qualitative detection of antibodies
to Hepatitis C Virus (HCV) in
human whole blood, serum or
plasma specimens



HASSLE-FREE

no special
equipment
needed



FAST

results in
20 minutes



EASY

sample preparation
protocol and visual
result interpretation



RELIABLE

test
performance



MINIMAL

sample
volume used

PRODUCT SPECIFICATION



Sensitivity
98.8%



Test time
15-20 minutes



Sample volume
50 µl (whole blood)
25 µl (serum /plasma)



Sample type
Human Whole Blood /
Serum / Plasma



Specificity
99.1%



Storage
2-30°C



Kit
25 tests



Shelf life
24 months

IVD

RESZON® HCV RAPID TEST



Whole Blood



Add 2 drops (approximately 50 µl) of whole blood sample into the sample well.

OR



Serum/Plasma



Add 1 drop (approximately 25 µl) of serum or plasma into the sample well.




Hold the buffer bottle vertically and add 1 drop of buffer (approximately 35 µl) into the sample well. Sample will start wicking up.



The results may be read from 15 to 20 minutes

! Allow the reaction to proceed for 15 minutes, read the results. Do not read results after 20 minutes.


RESULT INTERPRETATION

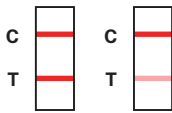

Positive



Negative

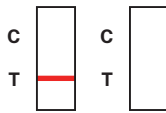

Invalid

REAGENTS AND MATERIALS SUPPLIED

- RESZON® HCV Rapid test cassette (individually sealed aluminium pouch with a desiccant)
- Assay buffer
- Optional: Disposable dropper, alcohol swab and safety lancet
- Instruction for Use (product insert)

ORDERING INFORMATION

RESZON® HCV RAPID TEST

Product Code	Description	Packing Size
RHC-RD0101-25	RESZON® HCV Rapid Test	25 tests / kit
RHC-RD0101A-25	RESZON® HCV Rapid Test (with disposable dropper)	25 tests / kit
RHC-RD0101B-25	RESZON® HCV Rapid Test (with disposable dropper, alcohol swab and safety lancet)	25 tests / kit

Contact us today to enquire: +6012 815 8689 | info@reszonics.com

Reszon Diagnostics International Sdn. Bhd. 201001033978 (917901-W)

A member of Hextar Healthcare Berhad

©2026 Reszon Diagnostics International Sdn. Bhd. All rights reserved. RESZON® is the registered trademark of Reszon Diagnostics International Sdn. Bhd. Doc No.: MM-01-04-37-EN; Rev. No.: 04; Effective Date: 2026.08.10
The content of this material is intended for Healthcare Professionals only, not for general public.

www.reszonics.com



Reszon Diagnostics

