

RESZON TYPHOID

RAPID TEST KIT

REF RTF-RD0101 / RTF-RD0202

MDA Reg. No.: IVDC85711391719

For professional *in vitro* diagnostic use only.

SCAN ME



Please scan the QR code for more information



RTF-RD0101

RESZON® Typhoid Rapid IgM

Detect specific IgM antibodies against a specific *Salmonella typhi* OMP antigen.



RTF-RD0202

RESZON® Typhoid Rapid IgG/IgM

Detect IgM and IgG antibodies against specific *Salmonella typhi* OMP antigen.



HASSLE-FREE

no any special equipment needed



FAST

results in 20 minutes



EASY

sample preparation protocol and visual result interpretation



RELIABLE

test performance



MINIMAL

sample volume used

RESZON TYPHOID

ENGLISH

RAPID TEST KIT

REF RTF-RD0101 / RTF-RD0202

MDA Reg. No.: IVDC85711391719

For professional *in vitro* diagnostic use only.




PRODUCT SPECIFICATION


	RTF-RD0101 RESZON® Typhoid Rapid IgM	RTF-RD0202 RESZON® Typhoid Rapid IgG/IgM
Sensitivity	85.9%	IgG: 94.44% IgM: 91.67%
Specificity	96.7%	IgG: 93.40% IgM: 95.24%
Storage condition	2-30°C	2-30°C
Specimen volume	35uL	45uL
Shelf life	24 months	24 months
Specimen	Human serum, plasma or whole blood	Human serum, plasma or whole blood
Assay time	15-20 mins (whole blood) 10-15 mins (serum/ plasma)	15-20 mins (serum/ plasma/ whole blood)

IVD


RESZON® TYPHOID RAPID IgG/IgM

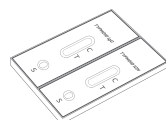

Label the test cassette with sample name

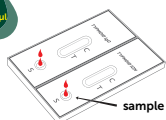

Add 45 µl whole blood/serum/plasma to each sample well.

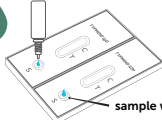


Add 1 drop of buffer to each sample well



Read results within 15-20 minutes

! Do not read the results after 30 minutes.


RESZON® TYPHOID RAPID IgM

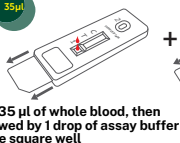


Add 35 µl of whole blood, then followed by 1 drop of assay buffer to the square well



Add 35 µl of serum/plasma to the square well (add 1 drop of assay buffer if sample does not wick up after 2-3 mins)

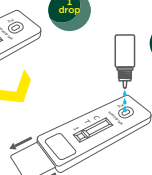

Let the sample flow until it reaches area C. Cassettes may be tapped gently on the table to facilitate sample flow

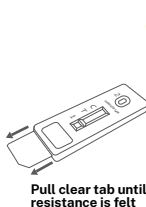

Add 3 drops of assay buffer to the oval well on top. Drop-by-drop and let it sip through the pad in between the drops



Read result after 15-20 mins for whole blood samples




Read result after 10-15 mins for serum/plasma samples


! Do not read the results after 30 minutes.

RESULT INTERPRETATION


NEGATIVE

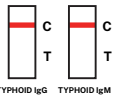

POSITIVE

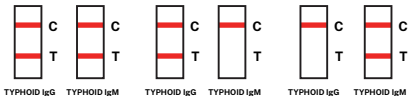

INVALID

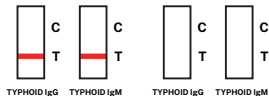

REAGENTS AND MATERIALS SUPPLIED

- RESZON® Typhoid Rapid IgG/IgM (Combo) cassette (individually sealed aluminium pouches)
- Assay buffer
- RESZON® Typhoid Rapid IgM cassette (individually sealed aluminium pouches)
- Instruction for Use (product insert)

ORDERING INFORMATION

RESZON® TYPHOID RAPID TEST KIT

Product Code	Description	Packing Size
RTF-RD0202	RESZON® Typhoid Rapid IgG/IgM (Combo) Test Kit	25 tests / kit
RTF-RD0101	RESZON® Typhoid Rapid IgM Test Kit	25 tests / kit

Contact us today to enquire: +6012 815 8689 | info@reszonics.com

Reszon Diagnostics International Sdn. Bhd. 201001033978 (917901-W)

A member of Hextar Healthcare Berhad

©2026 Reszon Diagnostics International Sdn. Bhd. All rights reserved. RESZON® is the registered trademark of Reszon Diagnostics International Sdn. Bhd. Doc No.: MM-01-04-03-EN; Rev. No.: 07; Effective Date: 2026.10.14
The content of this material is intended for Healthcare Professionals only, not for general public.

www.reszonics.com

Reszon
Diagnostics
INTERNATIONAL

Reszon Diagnostics

